

Bencox Bipolar Cup



Bencox Bipolar Cup




Optimal Fit

- A large array of component allows precise patient matching for optimal patient fit.
- Outer diameters from 38 to 60mm allow the surgeon to meet the unique challenges posed in revision and extended situations.
- The bipolar cup offers use with 22, 28mm femoral heads.

38 / 22	39 / 22	40 / 22	41 / 22	42 / 22	43 / 28	44 / 28	45 / 28	46 / 28	47 / 28	48 / 28	49 / 28
											
50 / 28	51 / 28	52 / 28	53 / 28	54 / 28	55 / 28	56 / 28	57 / 28	58 / 28	59 / 28	60 / 28	
											

Wide Range of Motion

- Enhanced 68° ROM resists the impingement to the UHMWPE locking ring.

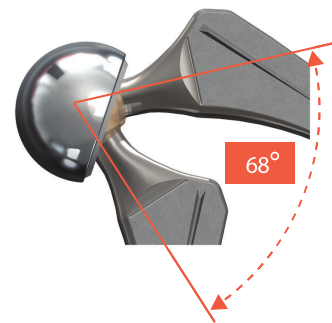

Proportional Eccentricity

- Multipolarity is achieved via multiple offsets within the system as cup sizes increase, the amount of eccentricity increases.
- Provides eccentricity that is proportional to the patient's anatomy.



38mm : 1mm

43mm : 1.5mm

50mm : 2.0mm

Simple Locking ring Mechanism



- Ring mechanism allows for easy assembly of the femoral head.



Enhanced Stability

- Superior pull-out strength resists the dislocation of the femoral head.

Testing Laboratory: ENDOLAB (No. 265.100430.10.1482)
Date: 26 May 2010

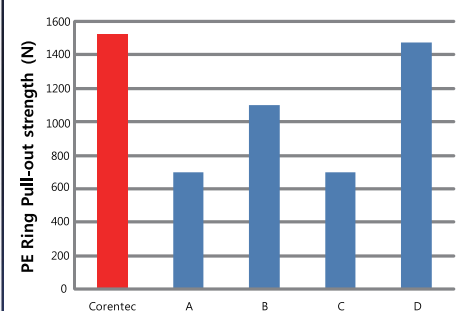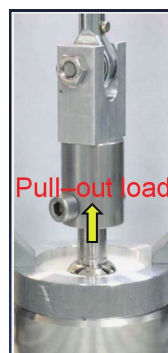

* Reference: Kummer FJ, Jaffe WL, Lesaka K, Perez F, Bull Hosp Jt Dis. 2005;63(1-2):5-8

Bencox Bipolar Cup •

Ordering Information

Part No.	Description
01.03.038	Bipolar Cup 38
01.03.039	Bipolar Cup 39
01.03.040	Bipolar Cup 40
01.03.041	Bipolar Cup 41
01.03.042	Bipolar Cup 42
01.03.043	Bipolar Cup 43
01.03.044	Bipolar Cup 44
01.03.045	Bipolar Cup 45
01.03.046	Bipolar Cup 46
01.03.047	Bipolar Cup 47
01.03.048	Bipolar Cup 48
01.03.049	Bipolar Cup 49
01.03.050	Bipolar Cup 50
01.03.051	Bipolar Cup 51
01.03.052	Bipolar Cup 52
01.03.053	Bipolar Cup 53
01.03.054	Bipolar Cup 54
01.03.055	Bipolar Cup 55
01.03.056	Bipolar Cup 56*
01.03.057	Bipolar Cup 57*
01.03.058	Bipolar Cup 58*
01.03.059	Bipolar Cup 59*
01.03.060	Bipolar Cup 60*

* Bipolar Cup 56, 57, 58, 59, 60
Custom Order (available upon request)

Manufacturing Plant 12, Yeongsanhong 1-gil, Ipjang-myeon, Seobuk-gu, Cheonan-si, Chungcheongnam-do, Korea
Tel : 82-41-585-7114 Fax : 82-41-585-7113

Seoul office 483, Gangnam-daero, Seocho-gu, Seoul, Korea Tel : 82-2-3445-5492~5 Fax : 82-2-3445-5467,5497

Corentec America 17305 Von Karman Ave. Suite 114, Irvine, CA 92614, USA Tel : 1-949-379-6227 Fax : 1-949-387-5716

EC Representative : Emergo Europe Molenstraat 15, 2513 BH, The Hague, The Netherlands Tel : 31-70-345-8570 Fax : 31-70-346-7299

<http://www.corentec.com>

CRSC-0100-M06(03)